
**Revisions**


Rev from	Rev to	ECO #
0504	0904	3050-04

**Notes:**


- BD Cat. Number 271040, 271041
- Blank (Sheet) Size : Length: 11"      Width: 25.5"  
 Number of Pages: 6      Number of Sheets: 1  
 Page Size: Length 11"      Width 8.5"      Final Folded Size: 2.75" x 4.25"
- Style (see illustrations below): # 4

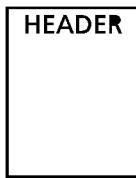
#1

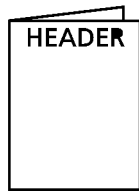
#2

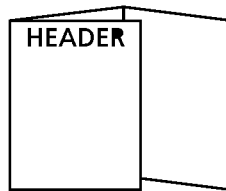
#3

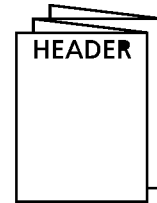
#4
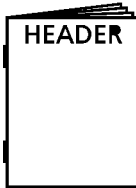
#5
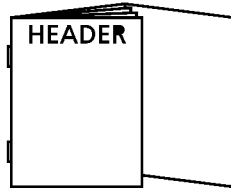
#6
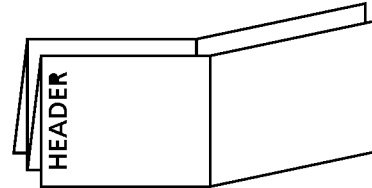
#7

- See Specification Control Number 8801151 for Material Information
- Ink Colors: Printed two sides  Yes     No  
 No. of Colors: 1      PMS# 2755 Blue
- Graphics are approved by Becton, Dickinson and Company. Supplier has the responsibility for using the most current approved revision level

VS Controlled by BD Caribe, LTD

Label Design	Date	COMPANY CONFIDENTIAL. THIS DOCUMENT IS THE PROPERTY OF BECTON, DICKINSON AND COMPANY AND IS NOT TO BE USED OUTSIDE THE COMPANY WITHOUT WRITTEN PERMISSION	 Becton, Dickinson and Company 7 Loveton Circle Sparks, MD 21152 USA	
Proofer	Date			
Checked By	Date			
Part Number:	8801151	Category and Description Package Insert, BBL GasPak Plus, Gas Generator Envelopes for Anaerobic	Sheet: 1 of 7 Scale: 1:1	<span style="font-size: 2em; font-weight: bold;">A</span>

# BD BBL™ GasPak™ Plus Anaerobic System Envelopes with Palladium Catalyst

CE 8801151  
2004/09

U.S. Patent No. 4,976,931

English: page 1 Italiano: pagina 4  
Français: page 2 Español: página 5  
Deutsch: Seite 3

See symbol glossary at end of insert. / Viz popis symbolů na konci příbalového letáku. / Se symbolglossaret i slutningen af indlægssedlen. / Zie lijst met symbolen aan het einde van de bijsluiter. / Vaadake sümbolite seletust infolehe lõpus. / Katso pakkauselosten lopussa olevaa kuvamerkkien sanastoa. / Voir le glossaire des symboles à la fin de la notice. / Siehe Symbol-Erklärungen am Ende der Packungsbeilage. / Δείτε το γλωσσάριο των συμβόλων στο τέλος του ένθετου. / A jelmagyarázat a használati utasítás végén található. / Vedere il glossario dei simboli alla fine del foglio illustrativo. / Ζr. informacino lapelio pabaigoje pateikiama simboli glosarijū. / Se i symbolforklaringen på slutten av produktvedlegget. / Zobacz objaśnienie symboli na końcu ulotki. / Consulte o glossário de símbolos no fim do folheto informativo. / Pozri slovník symbolov na konci letáka. / Consulte el glosario de símbolos al final del prospecto. / Se symbolförteckningen vid slutet av bipacksedeln.

Pokyny vám poskytne miestni zástupce spoločnosti BD. / Kontakt den lokale BD repræsentant for at få instruktioner. / Kasutusjuhiste suhtes kontakteeruge oma kohaliku BD esindajaga. / Επικοινωνήστε με τον τοπικό αντιπρόσωπο της BD για οδηγίες. / A használati utasítást kérje a BD helyi képviselőtől. / Naudojimo instrukcijų teiraukitės vietos BD įgaliotojo atstovo. / Aby uzyskać instrukcje użytkowania, skontaktuj się z lokalnym przedstawicielstwem BD. / Contacte o seu representante local da BD para obter instruções. / Instrukcie získate u miestneho zástupcu spoločnosti BD. / Kontakta lokal Becton Dickinson-representant för anvisningar.

## INTENDED USE

**GasPak™ Plus** is an improved disposable hydrogen and carbon dioxide generating envelope with catalyst sachet that produces an atmosphere in **BBL GasPak** jars suitable to support the primary isolation and cultivation of anaerobic bacteria.

## SUMMARY AND EXPLANATION

Prior to the mid-1960s, routine isolation and cultivation of anaerobes from clinical and other specimens was handicapped by the cumbersome methods necessary to achieve anaerobiosis in anaerobic jars.

In 1965, Brewer and Allgeier introduced a disposable hydrogen generator envelope<sup>1</sup>, which was modified to include carbon dioxide generation.

The performance of anaerobic bacteriology as a routine laboratory procedure was greatly facilitated when, in 1966, Brewer and Allgeier introduced a self-contained anaerobic system which eliminated the requirement for gas cylinders, vacuum pumps and other external equipment and, as an added safety feature, included an internal catalyst which did not require an external electrical connection.<sup>2</sup>

The **GasPak Plus** gas generator envelope differs from the **GasPak** envelope in that the catalyst pellets are in a sachet on the envelope. This method eliminates the need for catalyst pellet chambers which must be included in the jar lids of the **GasPak 100** and **GasPak 150** systems when standard **GasPak** envelopes are employed. The **GasPak Plus** envelope also features a unique easy open feature which allows easy access to the water channel of the envelope.

## PRINCIPLES OF THE PROCEDURE

The **GasPak Plus** envelope is activated by the addition of water, which passes through a series of channels to a filter paper wick. The wick feeds the water to the gas-generating tablets in the tablet chamber.

Hydrogen, generated from a sodium borohydride tablet following the addition of water, combines with the oxygen in the jar in the presence of the palladium catalyst to form water. Approximately 4 to 10% carbon dioxide, as measured by gas-liquid chromatography 1 h after activation, is generated from a sodium bicarbonate plus citric acid tablet. The carbon dioxide is provided to stimulate the growth of anaerobes that require or grow better in a CO<sub>2</sub>-enriched atmosphere.

The filter paper wick in the **GasPak Plus** envelope slows the introduction of water into the tablet chamber, allowing the lid to be placed on the **GasPak** jar before significant volumes of gases are released.

The performance characteristics of a properly operating **GasPak 100** system activated with a **GasPak** envelope have been described by Seip and Evans.<sup>3</sup>

Each **GasPak Plus** envelope has a catalyst sachet attached to the envelope, avoiding problems caused by inactivation of catalyst by hydrogen sulfide and other gases produced by some anaerobes.

## REAGENTS

Each **GasPak** disposable hydrogen + carbon dioxide generator envelope consists of:

- One Sodium Borohydride Tablet,
- One Sodium Bicarbonate Plus Citric Acid Tablet,
- One Filter Paper (piece),
- Palladium Catalyst Sachet.

## Warnings and Precautions

For *in vitro* Diagnostic Use.

May cause cancer by inhalation.

All usual precautions attendant to handling hydrogen gas should be observed when employing **GasPak Plus** hydrogen + carbon dioxide generator envelopes.

**Hydrogen is a flammable gas. A mixture of hydrogen gas with oxygen or air in a confined area will explode if ignited by a spark, flame or other source of ignition.**

The historical incidence rate of reported hydrogen gas ignition is approximately two (2) incidents per million envelopes.

## Catalyst

- Important: Do not allow the catalyst to come into contact with water. Careful addition of 10 mL of water is necessary.**
- THE APPROPRIATE **GasPak** SYSTEM RACK MUST BE UTILIZED AT ALL TIMES.
- GasPak Plus** envelopes with frays or cuts on the catalyst sachet covering should not be used.
- After use, the jar should be opened and allowed to aerate for approximately 15 sec prior to removing the rack and its contents so as to achieve atmosphere equilibration.

**BBL GasPak** hydrogen + carbon dioxide generator envelopes are intended for use only in properly maintained and operated **GasPak 100** and **GasPak 150** systems. Their use in other systems may yield misleading results and may even be hazardous to laboratory personnel. For satisfactory and safe use, carefully follow the instructions provided with each system. One envelope is to be used in each **GasPak 100** system and three envelopes are to be used in each **GasPak 150** system. Never mix **GasPak** and **GasPak Plus** envelopes in the same jar (i.e., envelopes with or without integral catalysts).

**ANY UNSATISFACTORY ENVELOPE MUST BE CUT OPEN AND EXAMINED FOR INTACT OR PARTIALLY INTACT TABLETS. ANY REMAINING TABLETS SHOULD BE PLACED IN A CONVENIENT SINK AND DISSOLVED IN FLOWING WATER. DO NOT DISCARD INTACT ENVELOPES.**

**Storage Instructions:** On receipt, store envelopes in a dry environment at 20 to 30°C. Avoid storage in close proximity to laboratory chemicals.

**GasPak Plus** envelopes are ready to use and should be activated according to instructions below and on the envelope. The expiration date is for unopened and intact envelopes.

Do not open until ready to use.

**Product Deterioration:** Do not use if envelopes or catalyst sachets show evidence of having been punctured, crushed, distorted, exposed to moisture, or exhibit other signs of deterioration.

## PROCEDURE

**Material Provided:** **GasPak Plus** Anaerobic System Envelopes with Palladium Catalyst (Cat. No. 271041; with indicators).

**Materials Required But Not Provided:**

- Culture media in Petri dishes, tubes, bottles or other suitable containers,
- Incubator (35°C),

- GasPak 100** system – Cat Nos. 260626, 260627 or **GasPak 150** system – Cat Nos. 260628, 260629,
- GasPak Dry** Anaerobic Indicator – Cat. No. 271051, unless provided with 271041
- Pipette, 10 mL, or syringe.

## Test Procedure:

**GasPak 100** and **GasPak 150** Anaerobic Systems

- Place inoculated plates or tubes (in **GasPak** tube holder) in the **GasPak 100** or **GasPak 150** rack. Peel open the corner of the envelope along the dotted line and place the envelope in the clip of the rack with the catalyst sachet toward the outside. One envelope is to be used in each **GasPak 100** system and three envelopes are to be used in each **GasPak 150** system. Do not crease, fold or crush the envelope.
- Open a **GasPak** dry anaerobic indicator and place in the small clip of the **GasPak** rack. Place the rack in the **GasPak** jar.
- Add 10 mL of tap, distilled or deionized water through the open/cut corner of the envelope. Use a pipette or syringe inserting only the tip into the open corner. Never push any object down into the envelope.
- Close the **GasPak** jar promptly after envelope is activated and finger tighten the lid clamp in the manner described in the insert directions for the system being used. Incubate at conditions appropriate for the organisms being cultured, but not above 42°C.
- After use, open the jar and allow it to aerate for approximately 15 sec prior to removing the rack with its contents. Do not handle the **GasPak** envelope until after the rack has been removed from the jar.

NOTE: Since the **GasPak Plus** envelope has an integral catalyst, the use of palladium catalyst pellets in the catalyst reaction chambers in the lids of **GasPak** jars is not required.

## User Quality Control:

- Each **GasPak** system should be tested periodically for its ability to provide adequate conditions for the growth of anaerobic bacteria.
- If visible condensation does not occur within 60 min after a **GasPak** envelope is activated, the reaction should be terminated by opening the jar. Check the condition of the lid and O-ring gasket as described below.
- The following quality control check list will help assure proper performance:
  - The system must contain fresh or rejuvenated catalyst if envelopes without integral catalyst are used.
  - The **GasPak** system operated with a **GasPak Plus** envelope, once activated, should be maintained upright at all times.
  - GasPak Plus** disposable hydrogen + carbon dioxide generator envelope – The envelope must be kept upright and peeled open along the dotted line at the upper right hand corner, as indicated on envelope. The gas-generating tablets must be present and not displaced or crushed, and the foil channels and reaction chambers must not be crushed or distorted and must be able to accept 10 mL of water. Improper positioning of envelope or use of a damaged envelope may lead to imbalances in reactivity. Damage to the channels, which regulate the reaction rate, may result in spattering and “boiling over” of the contents of the envelope and hence inadequate gas output. The envelope must be positioned in the jar with the catalyst facing outward.
  - Plastic lid and clamp – Inspect for cracks, chips, split O-ring grooves, and other irregularities. Most cracks (craze or stress lines) develop on prolonged use and forced or excessive tightening of clamps on lids. Early cracks may appear as fine black lines near the center of the lid or along the upper edges of the **GasPak 100** clamp near the screw. Cracked lids and clamps which exhibit these fine black lines should not be used.
  - O-ring gasket – Examine for cracks, bumps, dirt, and improper seating in the O-ring groove. Examine O-ring groove wall for splitting or cracking. The O-ring must be resilient, clean and snugly retained within the O-ring groove for satisfactory operation of the system. Replace as necessary.

Quality control requirements must be performed in accordance with applicable local, state and/or federal regulations or accreditation requirements and your laboratory's standard Quality Control procedures. It is recommended that the user refer to pertinent NCCLS guidance and CLIA regulations for appropriate Quality Control practices.

## RESULTS

Within 2 h of incubation at 35°C, the carbon dioxide concentration is greater than 4 but less than 10%. The **GasPak** dry methylene blue anaerobic indicator becomes decolorized within 6 h at 35°C.

## PERFORMANCE CHARACTERISTICS

In a study by Seip and Evans,<sup>3</sup> the oxygen concentration, carbon dioxide concentration and redox potential (Eh of plated media) were monitored using **BBL GasPak** Disposable Hydrogen + Carbon Dioxide Generator Envelopes. Sixty min after activation, the carbon dioxide concentration ranged from 4.6 to 6.5% and the oxygen concentration ranged from 0.2 to 0.7%. One hundred min after activation, the oxygen concentration had decreased to less than 0.2%.

## AVAILABILITY

- | Cat. No. | Description  |
|----------|--|
| 271040   | <b>BBL™ GasPak™ Plus</b> Anaerobic System Envelopes with Palladium Catalyst, Disposable H <sub>2</sub> + CO <sub>2</sub> Generator Envelopes, Pkg. of 10 envelopes   |
| 271041   | <b>BBL™ GasPak™ Plus</b> Anaerobic System Envelopes with Palladium Catalyst, Disposable H <sub>2</sub> + CO <sub>2</sub> Generator Envelopes, Pkg. of 100 envelopes (100 <b>GasPak</b> Dry Anaerobic Indicators also provided) |

## REFERENCES

- Brewer, J.H., and D.L. Allgeier. 1965. Disposable hydrogen generator. *Science* 147:1033-1034.
- Brewer, J.H., and D.L. Allgeier. 1966. Safe self-contained carbon dioxide-hydrogen anaerobic system. *Appl. Microbiol.* 14:985-988.
- Seip, W.F., and G.L. Evans. 1980. Atmospheric analysis and redox potentials of culture media in the **GasPak** system. *J. Clin. Microbiol.* 11:226-233.

# **BBL GasPak Plus Anaerobic System Envelopes with Palladium Catalyst**

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BENEX Limited