
Revisions

BALTSO0191 Version 10.0 Template 4
Inserts

Rev from	Rev to	JOB #
03	04	9056-17

NOTES:

- BD Catalog Number: 448876
- Blank (Sheet) Size: Length: 8.5" Width: 14.0"
- Number of Pages: 2 Number of Sheets: 1
- Page Size: Length: 8.5" Width: 14.0" Final Folded Size: No Fold
- Ink Colors: No. of Colors: 1 PMS#: Standard Black
- Printed two sides: Yes ☒ No ☐
- Style (see illustrations below): # 1

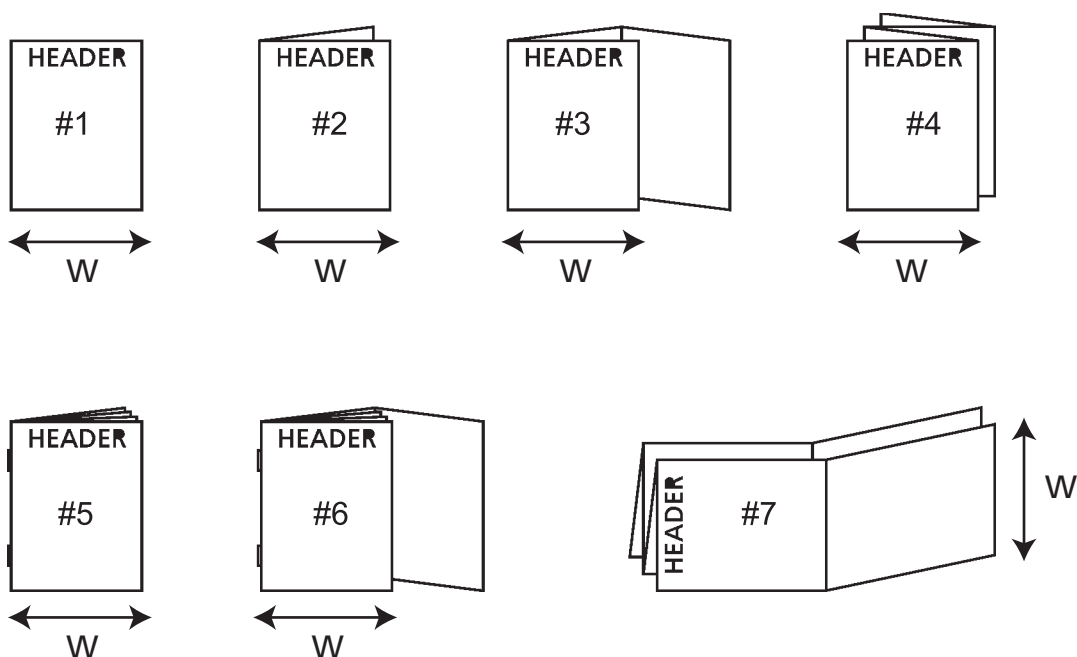- Vendor Printed ☐ Online/In House Printed ☒ Web ☐
- See specification control no. N/A for material information.
- Graphics are approved by Becton, Dickinson and Company. Supplier has the responsibility for using the most current approved revision level.

Label Design	REVISED BY By Tori Pagani at 1:17 pm, Mar 30, 2017	COMPANY CONFIDENTIAL. THIS DOCUMENT IS THE PROPERTY OF BECTON, DICKINSON AND COMPANY AND IS NOT TO BE USED OUTSIDE THE COMPANY WITHOUT WRITTEN PERMISSION.	 Becton, Dickinson and Company 7 Loveton Circle Sparks, MD 21152 USA
Proofer	PROOFING APPROVED BY By Natalie Morio at 12:18 pm, Apr 02, 2017		
Checked By	THIRD EYE BY By Terrence Means at 1:26 pm, Apr 10, 2017		
Part Number:	L010711	Category and Description Package Insert, Phoenix UNMIC-403 Panel	Sheet: 1 of 3 Scale: N/A A

INTENDED USE

The **BD Phoenix™** UNMIC panel is used for the susceptibility testing of most aerobic and facultative anaerobic gram-negative bacteria of human origin with selected antimicrobial agents (for a complete listing of taxa, refer to the **BD Phoenix** System User's Manual). This panel is only for use with the **BD Phoenix** Automated Microbiology System instrument.

PRINCIPLES OF THE PROCEDURE

Refer to the **BD Phoenix** System User's Manual.

PANEL

Susceptibility. The UNMIC panel contains the antimicrobial agents and concentrations in doubling dilutions found in Table 1.

STORAGE

Store at 15–25 °C. Do not use panel if the pouch is punctured or opened. Do not use the panel if desiccant is missing or if the desiccant pouch is torn.

Panels must be used within 2 h of being removed from the pouch.

Precautions:

For in vitro Diagnostic Use.

WARNINGS

Observe established precautions against microbiological hazards throughout all procedures. “Standard Precautions”^{1,2} and institutional guidelines must be followed in handling all items contaminated with specimens and microorganisms. Prior to discarding, sterilize specimen containers and other inoculated materials by autoclaving.

QUALITY CONTROL

Susceptibility. See expected results in Table 1.

NOTE: QC may be reported as less than or equal to the lowest, or greater than the highest concentration of the antimicrobic.

LIMITATIONS OF THE PROCEDURE

The clinical relevance of a specific antimicrobic and organism combination is determined by the activity of that antimicrobic against the organism and whether the antimicrobic is indicated for treatment of a disease state associated with that organism. The **BD Phoenix** System provides results for combinations, whether they are clinically relevant or not relevant. Overall, the results obtained using the **BD Phoenix** System compare favorably with the CLSI broth microdilution reference methodology;³ however, Table 2 contains clinically relevant combinations that did not meet the strict standards of BD during clinical trials and will either not be reported or an alternate method is recommended for confirmation of the result. Similarly, Table 3 contains combinations which are not clinically relevant and are either not reported or an alternate method is recommended for confirmation of the result.

PERFORMANCE CHARACTERISTICS

Refer to the **BD Phoenix** System User's Manual.

REFERENCES

- Clinical and Laboratory Standards Institute. 2005. Approved guideline M29-A3. Protection of laboratory workers from occupationally acquired infections, 3rd ed., CLSI, Wayne, Pa.
- U.S. Department of Labor, Occupational Safety and Health Administration, 29 CFR Part 1910.1030, 1991. Occupational exposure to bloodborne pathogens. Federal Register 56:64175-64182.
- Clinical and Laboratory Standards Institute. 2009. Approved standard M7-A8. Methods for dilution antimicrobial susceptibility tests for bacteria that grow aerobically, 8th ed., CLSI, Wayne, Pa.

Technical Information: In the United States contact BD Technical Service and Support at 1.800.638.8663 or www.bd.com.

ATCC is a trademark of the American Type Culture Collection.

© 2017 BD. BD, the BD Logo and BD Phoenix are trademarks of Becton, Dickinson and Company.

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Benex Limited