

# *IMPRESS*<sup>®</sup> system portfolio

Manage your instruments to help  
improve patient care

The *IMPRESS*<sup>®</sup> portfolio offers a full range of scalable, real-time and integrated instrument management solutions with extensive support. Using these solutions, you can manage your inventory, people and processes throughout the continuum of instrument management. The *IMPRESS* system can support your initiatives to:

- Enhance patient safety
- Reduce case delays
- Increase staff productivity
- Locate instruments and sets
- Limit immediate-use steam sterilization
- Control instrument spend
- Meet today's challenging industry standards



With more than 20 years of experience in instrument tracking, the V. Mueller<sup>®</sup> brand *IMPRESS* system is the name you can trust in the sterile processing department (SPD) and operating room (OR).




CareFusion

## Flexible, scalable instrument inventory management

Offering varying levels of functionality to meet your unique needs, the *IMPRESS* portfolio includes three easy-to-use, flexible solutions ranging from a basic system to an advanced, fully integrated system. All solutions feature the same look and feel, touch-screen technology, option to be web hosted or server deployed and multi-entity functionality. With the *IMPRESS* portfolio, you can easily upgrade to a system with more robust functionality as your needs evolve.

Features	<i>IMPRESS</i> Countsheets	<i>IMPRESS</i> Track	<i>IMPRESS</i> Connect
Countsheets management	•	•	•
350,000 instrument images	•	•	•
Catalog	•	•	•
Cross referencing	•	•	•
Reports	•	•	•
Inventory management		•	•
Process alerts		•	•
Washer in and out tracking		•	•
On-screen assembly		•	•
Critical item flagging		•	•
Repair management		•	•
Sterilization		•	•
Non-sterile warnings		•	•
Case cart management		•	•
Tracking to patient level		•	•
Warranty tracking		•	•
Advanced reporting		•	•
Connectivity/Interfaces			
– OR scheduler			
– Sterilizer			•
– Radio frequency identification (RFID)/Real-time locating systems (RTLS)			
Auto priority			•
Single instrument tracking			•
Loaner tracking			•
Immediate-use steam sterilization			•
OR reporting			•

### Track instruments throughout the continuum of instrument management

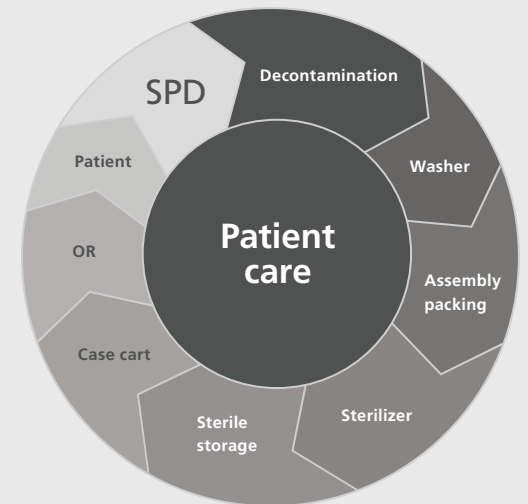

### *IMPRESS* system services and support

- Instrument inventory/validation
- Data entry
- 24/7 phone support
- Dedicated implementation managers
- Comprehensive V. Mueller instrumentation offering

## IMPRESS Countsheets system

IMPRESS Countsheets is the base-level offering of the IMPRESS portfolio. Through easy-to-use functionality, you can address your instrument management processes while managing information to drive operational efficiency.

### IMPRESS Countsheets features

- Automated countsheet management
- Countsheets standardization
- Countsheets access
- Reporting
- More than 350,000 instrument images
- V. Mueller catalog
- Cross reference tool

IMPRESS Countsheets enables you to manage, standardize and access countsheets

**MINOR BASIC TRAY**

Specialty:

Catalog Number	Manufacturer	Item Description	Std Qty	Actual Qty	OR 1st	OR Add	OR 2nd	OR Final
		<b>STRINGER</b>						
GL660	V Mueller	FOERSTER SPONGE FCPS STR SERR JAWS	1	0				
SU2936	V Mueller	NONPERF TOWEL CLAMP 5-1/4"	2	0				
SU2905	V Mueller	BACKHALS TOWEL FCPS 5-1/4"	6	0				
SU2702	V Mueller	MOSQUITO CVD SNAP	6	0				
SU2735	V Mueller	CRILE ART FCPS CVD 5-1/2"	10	0				
SU2700	V Mueller	PEAN ART FCPS 6-1/4" CVD	4	0				
SU2900	V Mueller	OCHSNER ART FCPS STR 6-1/4" 1X2 TEETH	4	0				
SU4055	V Mueller	ALLIS TISS FCPS 8" 5X8 TEETH	6	0				
SU5000	V Mueller	BABCOCK TISS FCPS 6-1/4"	2	0				
SU3900	V Mueller	LAHEY GOITER VUL FCP 8" STR 3X3 TEETH	4	0				
791851	V Mueller	VITAL METZ DISCT SCIS 5-3/4 IN CURVED REGULAR BSS	1	0				
MO1601	V Mueller	VITAL METZ FOR SCISS 7" CURVED*	1	0				
SU1814	V Mueller	VITAL MAYO 6-3/4" CVD DISCT SCISSORS	1	0				
SU1804	V Mueller	VITAL MAYO 6-3/4" STR DISCT SCISSORS	1	0				
792561	V Mueller	VITAL WEBSTER NEEDLE HOLDER 4-3/4" LT SMOOTH JAWS	2	0				
SU16060	V Mueller	VITAL MAYO-HEGAR NDL HLDR 6" *	2	0				
		<b>TOP OF TRAY</b>						
SU2303	V Mueller	DRESSING FORCEPS 6" *	2	0				
SU2333	V Mueller	TISSUE FCPS 1X2 TTH 6" *	2	0				
NL1405	V Mueller	ADSON FORCEPS 4-3/4" 1X2 TEETH*	1	0				
NL1410	V Mueller	ADSON FORCEP	1	0				
SU1405-001	V Mueller	SURGICAL HANDLE NO. 3 STAINLESS	1	0				
SU1407	V Mueller	SURGICAL HANDLE NO. 7 STAINLESS	1	0				
<b>Total Instrument Count</b>			<b>61</b>	<b>0</b>				

First Count \_\_\_\_\_  
 Circulator \_\_\_\_\_  
 Scrub \_\_\_\_\_  
 Assembler \_\_\_\_\_  
 Date / Time 08/31/2012 09:10 \_\_\_\_\_  
 Problem \_\_\_\_\_  
 Comment \_\_\_\_\_ Last Revised On: 08/20/2012 14:14 \_\_\_\_\_

Page 1 of 1

## IMPRESS Track system

In addition to providing all of the *IMPRESS* Countsheets capabilities, *IMPRESS* Track lets you easily track assets, manage sterilization and gain actionable intelligence for enhanced instrument management with robust functionality.

### IMPRESS Track features

- Multi-entity management
- Inventory management
- Process alerts
- Washer in and out tracking
- On-screen assembly
- Critical item flagging
- Repair management
- Sterilization management
- Non-sterile warnings
- Case cart management
- Patient-level tracking
- Warranty tracking
- 100+ reports
- Supervisor control center



## IMPRESS Track enables you to:

Proactively monitor and display SPD productivity for greater visibility.

The Supervisor Control Center interface displays a dashboard with several key sections:

- Inventory Status:** A table showing various instrument sets and their status.
 

Inventory Name	Last Reported	Last Scan	Status
MAJOR BASIC SET - 001	06/28/2012	06/29	71
MAJOR LARYNGOSCOPE 11447 - 001	06/28/2012	06/29	71
REPTD MICROSCOPE FOR 4437 - 001	06/28/2012	06/29	71
MAJOR LARYNGOSCOPE 11447E - 001	06/28/2012	06/29	71
DEPN-TOTAL-HEE (2208)-001	06/28/2012	06/29	71
- Inventory History:** A table showing historical data for the MAJOR BASIC SET - 001.
 

Inventory Name	Last Reported	Last Scan	Status
MAJOR BASIC SET - 001	06/28/2012	06/29	71
- Total Inventory Productivity:** A summary table showing overall performance metrics.
 

Area	Number Processed
Assemblies	2
Disassemblies	2
Scanners	2

View instrument and set images, indicate missing and damaged instruments, pull replacements from back wall inventory and pause the assembly process.

The Assembly interface provides a detailed view of a specific instrument set, including:

- Inventory Type:** MAJOR BASIC SET 1122 - 001
- Assembly Options:** View Backlog, Missing, Repair Scope, Replaced or Added, Replaced Scanned Set.
- Inventory Table:**

Item Name	Item Description	Std Qty	Min Qty
002704	PEARL ART FORCE CVD 8"	2	3
002708	PEARL ART FORCE 5-1/4" CVD (RESOLAP 5H430)	6	
002712	BARBAR ART FORCE CVD 5-1/4"	10	
002702	MOSGARD CVD SHAP	4	
002988	BACKHAUS TOWEL FORCE 5-1/4"	4	
002938	NOVA PERF TOWEL CLAMP 5-1/4"	2	
002989	COCHNER ART FORCE STR 5-9/16" 1/2 TEETH (RESOLAP 5H430)	6	
002984	COCHNER ART FORCE STR 5-9/16" 1/2 TEETH	2	
- Navigation:** Row Up, Row Down, Page Up, Page Down, Count Complete, Pause Count, Cancel Count, Dark Translucent.

Identify and flag items that are continuously in short supply for ongoing priority processing.

The Decontam interface displays a red alert box with the following text:

The following inventory items are scheduled for surgery. Please expedite them as a Priority.

Inventory Name:  
MAJOR BASIC SET - 001  
VASCULAR SET(GENERAL) 1147 - 001

Buttons: OK, Defect Inventory, Complete, Remove From List, Cancel Reading, Start Scanning, Exit.

Document sterilization and receive alerts when items are loaded into mismatched sterilizers and biological indicators need to be read.

The Sterilizer Station interface displays a red alert box with the following text:

**STERILIZER INVENTORY PARTS MISMATCH**  
This Item Cannot Be Sterilized in This Sterilizer:  
FLEXIBLE BRONCHOSCOPE ADULT 528 541 70437 - 001

AFTER THE ABOVE ITEM HAS BEEN REMOVED FROM THE LOAD SCAN IN YOUR STAFF CODE OR ENTER YOUR PIN NUMBER.

Barcode Staff Badge Or Enter Staff Pin

Buttons: 1, 2, 3, 4, 5, 6, 7, 8, 9, Cancel, 0, Clear, OK.

## IMPRESS Connect system

IMPRESS Connect offers the most advanced functionality in the IMPRESS portfolio, building on all of the IMPRESS Track capabilities. A fully interfaced and integrated open tracking solution that operates as a data and process repository, IMPRESS Connect enhances all sterile processing and distribution processes to help you optimize patient care.

### IMPRESS Connect features

- OR interface
- Sterilizer interface
- RFID/RTLS connectivity
- Case cart management
- Loaner tracking
- Single instrument tracking
- Robust reporting





To learn more about the *IMPRESS* portfolio and solutions that meet your unique needs, request a **workflow/process assessment** today. Contact your V. Mueller sales representative, call Customer Service at **800.323.9088** or visit [carefusion.com/vmueller](http://carefusion.com/vmueller).

CareFusion  
Waukegan, IL

[carefusion.com](http://carefusion.com)

