



Indispensable to  
human health

# **DrySlide™**


*The Dry Alternative  
In Rapid Testing*

Now there's a dry alternative for rapid testing with microbiological biochemical reagents. Why use awkward disks, strips, filter papers, messy liquids or vials of reagents? Now you have the option of the disposable, convenient and cost-effective **DrySlide™** family of products. This innovative system employs a dry reagent cassette for spot biochemical analysis of bacteria. And because it's divided into four large segments, you can do multiple testing on a single slide. There's no spillage or breakage to worry about. Plus, the device is totally self-contained. And facilities that perform these tests frequently will find that **DrySlide™** is more cost-effective than other methods. The **DrySlide™** line is available in the following recommended methodologies:



### DrySlide™ Indole

A disposable slide containing a reagent cassette impregnated with 5% *p*-dimethylaminobenzaldehyde (DMABA) for the rapid "spot" identification of indole producing bacteria from primary culture. Test results in 30 seconds.

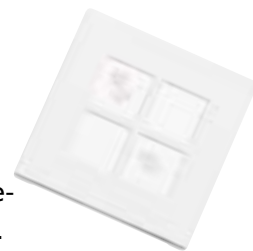

### DrySlide™ Oxidase

A disposable slide containing film coated with oxidase reagent (N,N,N',N'-tetramethyl-*p*-phenylenediamine dihydrochloride). Test results in 20 seconds.



### DrySlide™ Nitrocefin

A disposable slide containing a reagent cassette impregnated with nitrocefin for the rapid qualitative detection of beta-lactamase producing bacteria. Test results in 60 minutes or less.

### DrySlide™ PYR Kit

**DrySlide™** PYR is a disposable slide format for the presumptive identification of group A streptococci and enterococci. Contains the reagent L-pyrrolidonyl- $\beta$ -naphthylamide (PYR) which is hydrolyzed by bacteria possessing the enzyme pyrrolidonyl amino peptidase. The addition of the **DrySlide™** PYR Color Developer, which is provided in convenient crush ampules, results in a vivid color reaction in a 3 minute test.



**BD Biosciences**  
7 Loveton Circle  
Sparks, MD 21152  
tel: 410.316.4000  
www.bd.com

Product	Description	Original Cat #	New Codified No.
DrySlide Indole	75 slides	4331748	231748
DrySlide Oxidase	75 slides	4331746	231746
DrySlide Nitrocefin	25 slides	4331749	231749
DrySlide PYR Kit	15 slides - 15 ampules	4331747	231747