

## Trust the BD Vacutainer® Brand— A System of products you can rely on

At BD, we pledge to continually support the healthcare professional's need for quality educational tools and training assistance with practical products and professional development materials.



### Ordering Information

#### BD Vacutainer® Push Button Blood Collection Sets with Pre-Attached Holder

Reference Number	Needle Gauge	Wing Color	Tubing Length	Configuration	Packaging
367352	21	21G	12"	with holder	20/Box 100/Case
368656	23	23G	12"	with holder	20/Box 100/Case
368659	25	25G	12"	with holder	20/Box 100/Case

#### BD Vacutainer® Push Button Blood Collection Sets

367338	21	21G	7"	with luer	50/Box 200/Case
367344	21	21G	12"	with luer	50/Box 200/Case
367326	21	21G	12"	without luer	50/Box 200/Case
367336	23	23G	7"	with luer	50/Box 200/Case
367342	23	23G	12"	with luer	50/Box 200/Case
367324	23	23G	12"	without luer	50/Box 200/Case
367341	25	25G	12"	with luer	50/Box 200/Case
367323	25	25G	12"	without luer	50/Box 200/Case

BD Technical Services: **1.800.631.0174** (US only)  
 BD Customer Service: **1.888.237.2762** (US only)  
 BD Customer Service: **1.866.979.9408** (Canada only)  
 or visit [www.bd.com/vacutainer](http://www.bd.com/vacutainer)



**BD Diagnostics**  
 Preanalytical Systems  
 1 Becton Drive  
 Franklin Lakes, NJ 07417  
[www.bd.com/vacutainer](http://www.bd.com/vacutainer)

<sup>1</sup> Mary Hotaling: "A Retractable Winged Steel (Butterfly) Needle Performance Improvement Project." *Joint Commission Journal on Quality and Patient Safety*, Volume 35, Number 2, pages 100-105, 2009.  
<sup>2</sup> OSHA Safety and Health Information Bulletin (SHIB 10-15-03). Accessed at: [www.osha.gov](http://www.osha.gov).  
<sup>3</sup> National Phlebotomy Association Safety Alert. Accessed at: [www.nationalphlebotomy.org](http://www.nationalphlebotomy.org)



► **Safety**  
 from End to End



## INTRODUCING

**BD Vacutainer®**  
 Push Button Blood Collection  
 Set with Pre-Attached Holder



Helping all people  
 live healthy lives

# Safety from End to End

## Needle Safety

The **BD Vacutainer® Push Button Blood Collection Set with Pre-Attached Holder** helps guard against costly, potentially debilitating needlesticks, allowing caregivers to focus on superior patient care.

### Safety with the push of a button


- Provides split-second protection against needlestick injury
- Clinically demonstrated to reduce accidental needlestick injuries by **88%** with zero injuries in the final 21 months of the study<sup>1</sup>
- Ease of use reduces training requirements and is designed to improve activation compliance
- Mitigates significant costs associated with needlestick infections

### One-handed safety activation

- Allows for easy activation of the safety mechanism while attending to the patient/venipuncture site
- Intuitive design makes blood collection safer without compromising patient care

### In-vein needle retraction for a positive patient experience

- Reduces the risk of professional or patient exposure to a contaminated needle
- Delivers easy activation while minimizing patient discomfort
- Telling patients that the audible “click” signifies the procedure is complete may improve patient experience

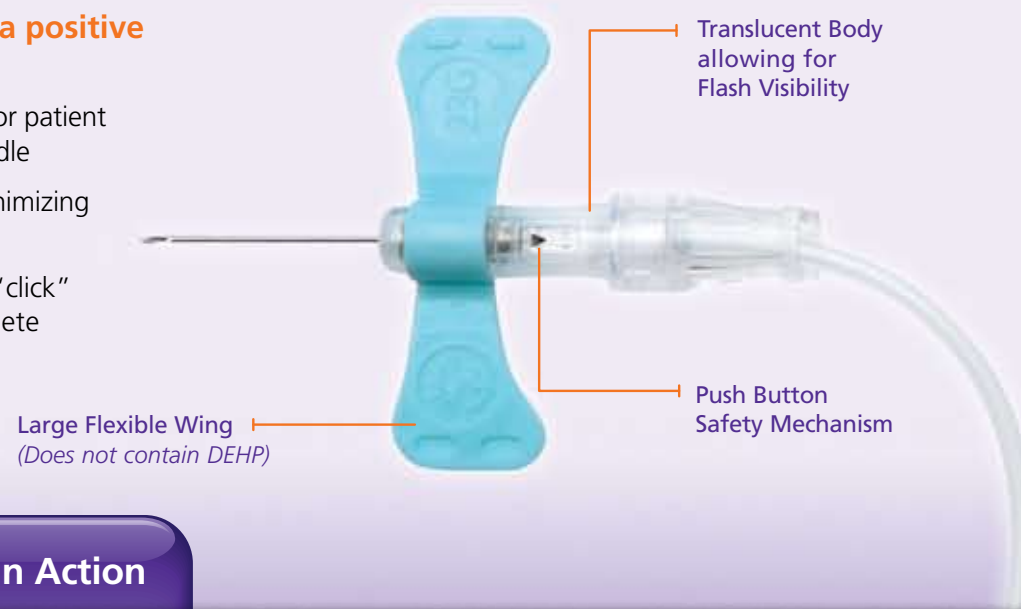

## Holder Safety


The **BD Vacutainer® Push Button Blood Collection Set with Pre-Attached Holder**

- ▶ **Safety**  
Pre-attached holders minimize the risk of accidental needlestick injuries
- ▶ **Compliance**  
Pre-attached holders are compliant with OSHA's single-use mandate<sup>2</sup>
- ▶ **Efficiency**  
Pre-attached holders eliminate steps in the phlebotomy process

### Now Available with Convenient Pre-Attached Holder

According to a recent National Phlebotomy Association study, 99% of the blood collection tube holders evaluated had been contaminated with blood — putting both sample quality and worker safety at risk.<sup>3</sup>

In 2003, an OSHA Safety and Health Bulletin prohibited the re-use of tube holders. **BD Vacutainer® Push Button Blood Collection Sets** are equipped with pre-attached disposable holders to protect caregivers from back-end needlesticks, offering an even higher level of safety.<sup>2</sup>

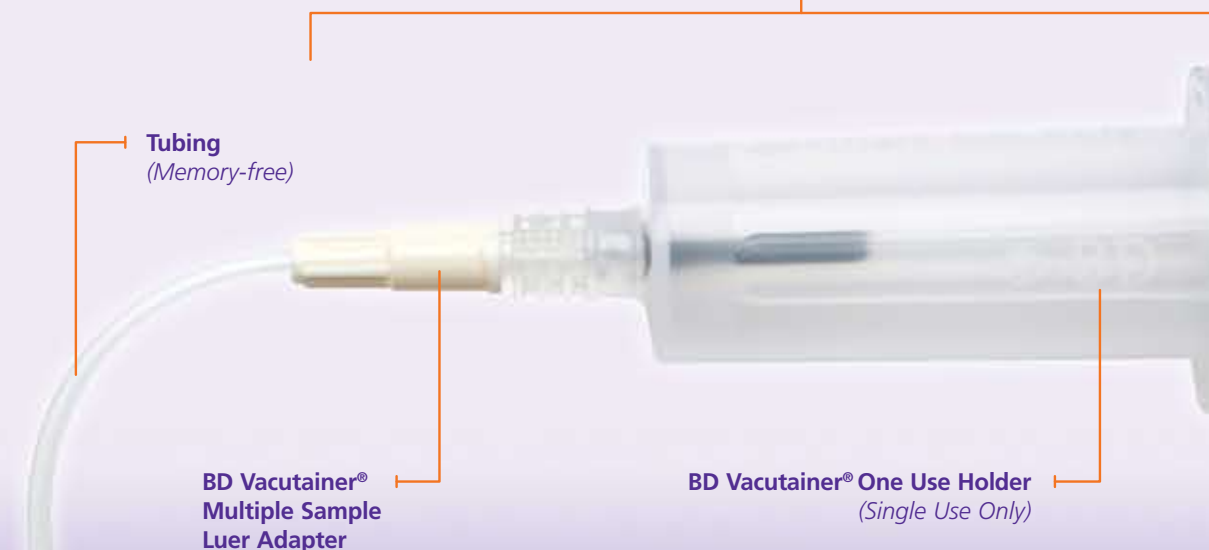

## Safety & Efficiency in Action



Memory-free tubing:  
Tubing does not coil or curl



Before activation



After activation



Easy & flexible: Choose the technique that works for you



Flash visualization confirms venous access



Easy in-vein activation without additional patient discomfort



Needle retracts from the vein and locks into place, offering immediate protection against needlestick injury

View the **BD Vacutainer® Push Button Blood Collection Set** in action and learn more about this innovative advance in safety and satisfaction:  
[www.bd.com/vacutainer/producttraining](http://www.bd.com/vacutainer/producttraining)