


Channel error

BD Alaris™ System

What does channel error mean?

The BD Alaris™ modules run self-checking programs prior to and during operation. A **CHANNEL ERROR** message means the device has discovered an error either in the affected channel hardware or software. **Operation on the affected channel stops; other infusing channels will not be affected.**

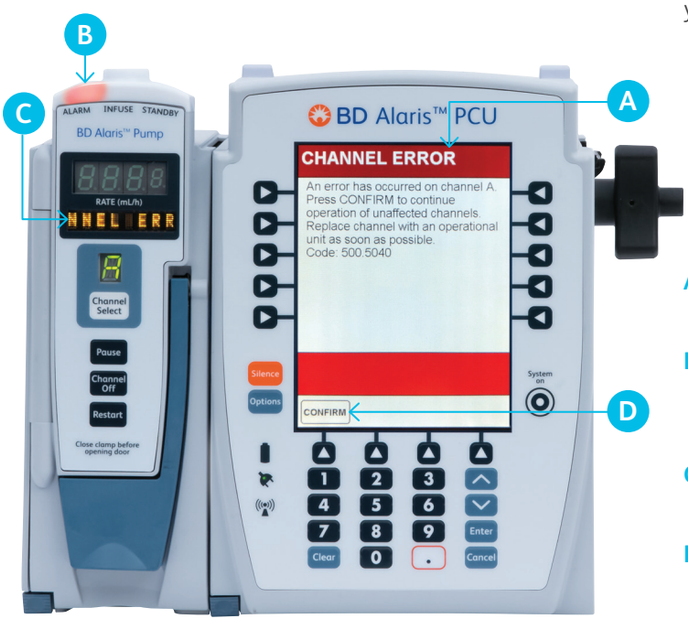

What should I do?

To silence the alarm and continue unaffected channel operation, press the **CONFIRM** soft key. Obtain a new module. Operation of the affected channel stops. It may be necessary to utilize the restore function on unaffected channels that may have been attached and reattached to the BD Alaris™ PC unit during the alarm state. Return the affected module to your facility's biomed department.

A. The system will display a **CHANNEL ERROR** message on the main display and provide the user with instructions.

B. The flashing red light on the module and the audible alarm tones will indicate that the unit is in an alarm state. The affected channel will stop operation.

C. The **CHANNEL ERROR** message will scroll on the affected module's display.

D. The **CONFIRM** soft key will appear at the bottom left of the main display.

For product support, contact Customer Advocacy at **888.812.3266** or **customerfeedback@bd.com**

For technical support, contact Instrument Technical Support at **888.812.3648.**

For product orders, contact Customer Order Management at **800.482.4822.**

⚠ See reverse side for applicable warnings and cautions.

For complete instructions, refer to the BD Alaris™ System User Manual at **bd.com**

Channel error

⚠ Warnings and cautions

WARNING: The BD Alaris™ System performs a self check during power up. The PC unit should beep, no errors should occur, and if a module is connected, all LED segments should flash. If the BD Alaris™ System fails the self check, remove the failing PC unit or module from use.

BD, San Diego, CA, 92130, U.S.

bd.com

BD, the BD Logo and Alaris are trademarks of Becton, Dickinson and Company or its affiliates.
© 2019 BD. All rights reserved. MC190002458 (0819/500) BD-11199